

# PT01-ProteanFect™ Transfection Kit: Instructions for Use

## Product Overview

ProteanFect™ Transfection Kit offers a non-viral, non-electroporation, and non-liposomal transfection system utilizing engineered mammalian proteins. This innovative design achieves high transfection efficiency while maintaining a superior safety profile. Specifically developed for hard-to-transfect cell lines, the kit ensures robust performance across a broad range of cell types (refer to Table 4). Additionally, it excels in delivering high molecular weight nucleic acids and co-transfecting multiple nucleic acids.

## Component Description

The kit is shipped on dry ice. Once received, store the components as indicated below. The kit includes positive control samples with EGFP-encoding mRNA and pDNA (plasmid DNA) to verify transfection efficiency.

**Table 1 Storage Conditions for the Components**

Component	Storage
Reagent A	2-8°C
Reagent B	-20°C
EGFP mRNA (1 µg/µL)	-80°C
GFP pDNA (0.5 µg/µL)	-20°C

**Note:** Avoid repeated freeze-thaw cycles of Reagent B, EGFP mRNA, and EGFP pDNA.

## Pre-Experimental Preparation

**Cell Condition:** Ensure cells are in optimal physiological condition on the day of transfection, with >90% viability.

**Reagent:** Allow Reagents A and B to reach room temperature. Briefly mix each reagent by inverting or vortexing prior to use.

**Medium:** Use Opti-MEM as the recommended medium. Serum-free RPMI 1640 or DMEM can be used as alternatives. Pre-warm the medium to 37°C or room temperature before use.

## PT01-ProteanFect™ Transfection Kit: Instructions for Use

### Transfection Procedure

**Table 2 Transfection Protocol for mRNA per Well of a 96-Well Plate**

Steps	Instructions for Cell Lines
<b>1. Transfection Complex Preparation <sup>a</sup></b>	
1.1 Mix Reagent A with mRNA	Mix 0.5 µg of mRNA with 40 µL of Reagent A. <b>Note:</b> Invert the Reagent A briefly before use to ensure uniformity.
1.2 Add Reagent B	Add 1.4 µL of Reagent B to the mixture. Mix gently by pipetting up and down 20-30 times or vortexing for 10 seconds.
<b>2. Cell Preparation</b>	
2.1 Suspension cells	Harvest the cells by centrifugation at 300 g for 5 minutes. Discard the supernatant and wash cells once with Opti-MEM. Resuspend cells with Opti-MEM and adjust concentration to 5×10 <sup>6</sup> - 1×10 <sup>7</sup> cells/mL. <b>Note:</b> Avoid including FBS in the transfection medium.
2.2 Adherent cells	Maintain 50%-80% cell confluence. Remove medium, wash cells once with Opti-MEM, then add 20 µL of Opti-MEM. <b>Note:</b> Avoid including FBS in the transfection medium. <b>Optional:</b> Harvest cells by trypsinization, then resuspend them in Opti-MEM at a concentration of 5×10 <sup>6</sup> - 1×10 <sup>7</sup> cells/mL for subsequent transfection.
<b>3. Transfection</b>	
3.1 Mix transfection complex with cells	For suspension cells, mix 40 µL of transfection complex with 20 µL of cell suspension and gently pipet up and down 2-3 times. For adherent cells, apply directly to the cells.
3.2 Incubation	Incubate the cells with the transfection complex for 45-60 minutes in a cell culture incubator.
3.3 Termination	Terminate the reaction by adding ≥200 µL of culture medium (10X cell suspension), centrifuge at 300 g for 5 minutes, and discard the supernatant. For adherent cells, replace the transfection mixture with ≥200 µL of culture medium (10X cell suspension). <b>Note:</b> The cell pellet may adhere to the tube walls. Gently remove the supernatant to minimize cell loss.
3.4 Post-transfection culture	Incubate the transfected cells in culture medium and assess transfection efficiency after 5 to 48 hours, or at an appropriate time.

FBS, Fetal bovine serum. **a.** The transfection complex may become slightly viscous during preparation. It can be directly added to cells once prepared without incubation. For optimal results, use the complex within 30 minutes.

## PT01-ProteanFect™ Transfection Kit: Instructions for Use

**Table 3 Transfection Guidelines for Different Culture Formats**

Components	Culture Vessels <sup>a</sup>	Cell Lines		
<b>Reagent A</b>	96-well	40 µL		
	48-well	80 µL		
	24-well	200 µL		
	12-well	600 µL		
	6-well	800 µL		
<b>Nucleic Acids <sup>b</sup></b>		<b>DNA</b>	<b>mRNA</b>	<b>siRNA</b>
	96-well	0.4 µg	0.5 µg	40 pmol
	48-well	0.8 µg	1 µg	80 pmol
	24-well	2 µg	2.5 µg	200 pmol
	12-well	6 µg	7.5 µg	600 pmol
<b>Reagent B</b>	96-well	1 µL	1.4 µL	1.4 µL
	48-well	2 µL	2.8 µL	2.8 µL
	24-well	5 µL	7 µL	7 µL
	12-well	15 µL	21 µL	21 µL
	6-well	20 µL	28 µL	28 µL
<b>Recommended Cell Number (Opti-MEM) <sup>c</sup></b>	96-well	$1 \times 10^5 \sim 2 \times 10^5$ (20 µL)		
	48-well	$2 \times 10^5 \sim 4 \times 10^5$ (40 µL)		
	24-well	$5 \times 10^5 \sim 1 \times 10^6$ (100 µL)		
	12-well	$1.5 \times 10^6 \sim 3 \times 10^6$ (300 µL)		
	6-well	$2 \times 10^6 \sim 4 \times 10^6$ (400 µL)		

**a.** For large-scale transfections, such as in 48-well plates or larger formats, it is recommended to use centrifuge tubes for the transfection process. **b.** When co-transfecting multiple nucleic acids, please ensure the total amount of nucleic acids added matches the recommended quantities for each plate format, as outlined in Table 3. **c.** The recommended cell number is primarily for suspension cells. For adherent cells, please adjust the cell number based on confluency.

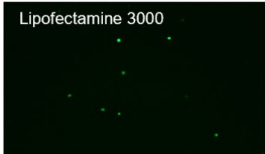
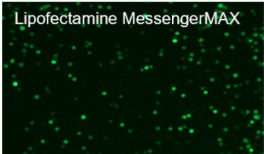
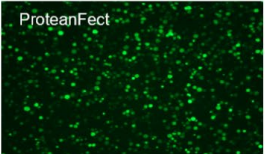
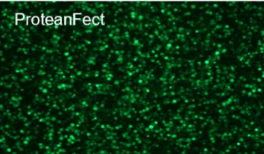
## PT01-ProteanFect™ Transfection Kit: Instructions for Use

**Table 4 Cell Lines Successfully Transfected Using ProteanFect™ Transfection Kit**

Cell Lines	Tested Nucleic Acid Types for Transfection
Jurkat T (human T lymphoblastic leukemia cells)	pDNA, mRNA, siRNA
LX-2 (human hepatic stellate cells)	pDNA, mRNA, siRNA
HepG2 (human liver tumor cells)	pDNA, mRNA, siRNA
THP-1 (human acute monocytic leukemia cell line)	mRNA, siRNA
Raji (human Burkitt's lymphoma cells)	mRNA, siRNA
K562 (human chronic myeloid leukemia cells)	pDNA, mRNA, siRNA
MOLT-16 (human T lymphoblastic leukemia cells)	mRNA, siRNA
SH-SY5Y (human neuroblastoma cells)	pDNA, mRNA, siRNA
U2OS (human osteosarcoma cells)	pDNA, mRNA, siRNA
U937 (human lymphoma cell line)	mRNA, siRNA
HFF (human foreskin fibroblasts)	pDNA, mRNA, siRNA
HEK-293 (human embryonic kidney cell line)	pDNA, mRNA, siRNA
MC38 (mouse colon cancer cells)	mRNA, siRNA
RAW264.7 (mouse mononuclear macrophage leukemia cells)	mRNA, siRNA
LLC (mouse Lewis lung cancer cells)	pDNA, mRNA, siRNA
C2C12 (mouse myoblasts)	pDNA, mRNA, siRNA
COS7 (African green monkey kidney fibroblast-like cells)	pDNA, mRNA, siRNA

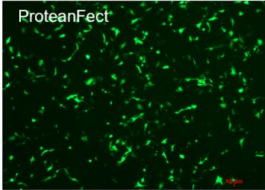
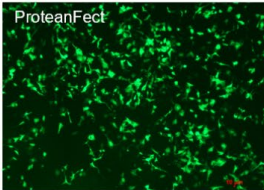
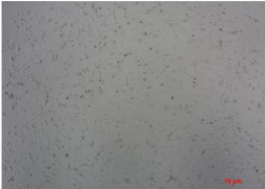

## Supporting Data

**Case 1: Successful Transfection of Jurkat T Cells**

EGFP DNA	EGFP mRNA
	
	

- **Cell Type:** Jurkat T Cells (suspension cell)
- **Transfected Nucleic Acid Type:** EGFP pDNA or EGFP mRNA
- **Detection Time:** 48 hours post-pDNA transfection; 24 hours post-mRNA transfection

**Case 2: Successful Transfection of LX-2 Cells**

1E4-1h	2.5E4-1h
	
	

- **Cell Type:** LX-2 Cells (adherent cell)
- **Cell Density:** The cells were plated at a density of  $1 \times 10^4$  cells per well and  $2.5 \times 10^4$  cells per well one day prior to transfection
- **Transfected Nucleic Acid Type:** EGFP pDNA
- **Detection Time:** 48 hours post-pDNA transfection

## Frequently Asked Questions (FAQs) and Troubleshooting Guide

### 1. Low Transfection Efficiency

#### 1.1 Optimize Transfection Parameters

Optimize transfection parameters for each cell type. **Extend incubation time:** Adjust the incubation time of the transfection complex with cells. The maximum incubation time is 2 hours for cell lines, and 30 minutes for primary cells. **Increase ProteanFect transfection complex:** Consider increasing the amount of transfection complex to improve transfection efficiency.

#### 1.2 Improve Plasmid DNA quality

Transfection efficiency is highly dependent on the quality of the plasmid DNA. **Endotoxin-free:** Use an endotoxin-free kit to prepare plasmid DNA. **Optimal concentration:** Ensure an OD<sub>260/280</sub> ratio between 1.7 and 1.9, and dilute DNA to 0.5–2 µg/µL using nuclease-free water.

#### 1.3 Improve Cell Condition

For cell lines, transfect cells with >90% viability, confirmed by trypan blue exclusion. Avoid using cells beyond 15 passages, and allow 2-3 passages for recently thawed cells to stabilize before transfection.

#### 1.4 Use Positive Control

We recommend using a 96-well plate format to optimize transfection conditions for a specific cell type, with EGFP mRNA as the positive control.

### 2. Low Cell Viability

Transfected cells may exhibit transient changes in behavior, but typically, viability should be restored by the second day post-transfection.

### 3. Lack of Pellet Post-Centrifugation

In 96-well formats, it is common for the pellet to be less distinct and may adhere to the tube walls. Gently pipetting can help minimize cell loss.

**Contact Information:** For further questions, please contact us at: [proteanflect@nanoportlabio.com](mailto:proteanflect@nanoportlabio.com).